

Accompanying the Recovery Journey

Create a Safer Environment

INFANT
RESUSCITATOR
www.procaremedfurn.com

KEY FEATURE

- **Reliable Performance:**
- No Power or Battery requirement, our Infant Resuscitator offers a silent operation. It delivers a consistent Peak Inspiratory pressure (PIP) and Positive End Expiratory Pressure (PEEP). Ensuring the proper and efficient ventilation for Infant
- **Connectivity:**
- The patient T-Piece can connect easily to neonatal masks or endotracheal tubes directly. With different sizes of face masks catering to vrying baby sizes, our resuscitator accommodates infants of all dimensions, ensuring a snug and secure fit.
- **Safety Features:**
- Our Infant Resuscitator is equipped with an inspiratory high pressure safety cut-off relief setting, and it prioritizes the safety of the newborn babies at all times. The adjustable PIP and PEEP setting allow doctors to tailor ventilation to individual needs, ensuring care with every use.
- **Efficient Operation:**
- It works on low flow from 5 to 15 liters per minute. It maximizes efficiency while minimizing resource Consumption. The inclusion of a medical-grade manometer enables real time monitoring of delivered PIP and PEEP, empowering healthcare providers with accurate data for informed decision making.
- try our other [Neonatal Care Machines](#)
-

Grab
Handle



Adjustable
Knob




Pipe



CE

ProCare

Address: Gate no 197, Dehu Alandi road,
Jadhavwadi, Chikali Pune,
Maharashtra-411062
Mob no: +91 8377832832
Web Site : www.procaremedfurn.com
GST NO: 27ABCFP8176F1ZY

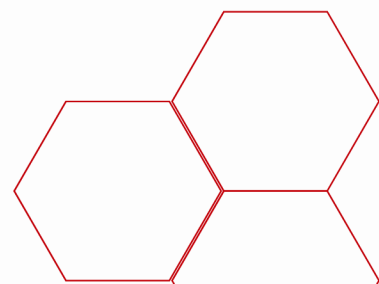

Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

ProCare