

# Accompanying the Recovery Journey

Create a Safer Environment

SUCTION MACHINE  
WITH GLASS JAR  
SKU: PRO192

[www.procaredmedfurn.com](http://www.procaredmedfurn.com)

## KEY FEATURE

### STRONG SUCTION:

THIS SUCTION MACHINE ALLOWS YOU TO CONTROL THE SUCTION, THE MACHINE COMES WITH A POWERFUL  $700 \pm 10$  MM HG SUCTION CAPACITY AND IT IS ADJUSTABLE TO 20-30 LITERS PER MINUTE

### EXTRAORDINARILY QUIET OPERATION:

THE MACHINE IS DESIGNED FOR SILENT AND EFFICIENT WORK, WHICH WILL HELP MAINTAIN THE SILENCE WITHIN THE WARDS.

### POWER:

THIS SUCTION MACHINE NEEDS 230V, 50 HZ,  $4 \pm 0.5$  AMPS, AND 200 WATTS TO PROVIDE CONSISTENT POWER FOR A SMOOTH WORKING

### ACCURATE MEASUREMENT:

THE SUCTION MACHINE WITH GLASS JARS IS EQUIPPED WITH A 2-INCH BCPM VACUUM GAUGE WHICH IS CALIBRATED BETWEEN 0-760 MM HG, AND IT HELPS TO MANAGE THE ACCURATE SUCTION LEVEL AND ADJUST IT ACCORDING TO THE PROCEDURES.

### HIGH-QUALITY PUMP:

TO ANY SUCTION MACHINE, ITS MOTOR AND PUMPS ARE IMPORTANT FOR ACCURATE AND EFFICIENT WORK, THIS MACHINE HAS AN ELECTRIC MOTOR OF 1/4 OR 1/2 HORSEPOWER AND AN OIL-IMMERSED PISTON PUMP FOR LONG OPERATIONS.

### STURDY GLASS JARS:

THIS SUCTION MACHINE COMES WITH TWO GLASS JARS, THE CAPACITY OF EACH JAR IS 1.5 LITER. THE RUBBER COVERS ARE PROVIDED FOR THE JARS THAT SHUT THEIR OWN FOR SAFE AND SECURE CONTAINMENT OF LIQUID

### ROBUST HOUSING:

THE VACUUM SUCTION MACHINE WITH A GLASS JAR IS SECURED IN A POWDER-COATED MILD STEEL CABINET EQUIPPED WITH A GRAB HANDLE, AND RUBBER BUSHING FOR BETTER STABILITY.

### SMALL AND EASY TO TRANSPORT:

THIS SUCTION MACHINE IS LIGHTWEIGHT AND EASY TO CARRY, THE SIZES ARE 19 X 12 X 29 INCHES, AND THE WEIGHT OF THE MACHINE IS ONLY 19 KG.

Jar



Grab  
Handle




Wheel



CE

# ProCare

Address: Gate no 197, Dehu Alandi road,  
Jadhavwadi, Chikali Pune,  
Maharashtra-411062  
Mob no: +91 8377832832  
Web Site : [www.procaredmedfurn.com](http://www.procaredmedfurn.com)  
GST NO: 27ABCFP8176F1ZY

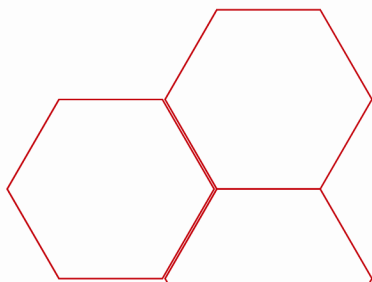

Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

ProCare