

Accompanying the Recovery Journey

Create a Safer Environment

YUWELL 7E-G1 SUCTION MACHINE SKU: PRO199 www.procaredmedfurn.com

The required input power source for 7E-

G1: **DC 12 V / 90 VA / AC 220 V @ 50 Hz**

Comes with a **12 V lithium-ion battery**

Pumping Rate of the Machine: **25 L/min**

80 kPa is the maximum negative pressure.

Noise level during operation: less than **50 dB**

Capacity of Jar: **1000 mL**

Operating Duration: **60 minutes** when
using a battery

KEY FEATURE

DUAL POWER SOURCE:

THE 7E-G1 PORTABLE SUCTION MACHINE IS SUITABLE TO USE IN ANY CONDITION YOU CAN USE THE MACHINE ANYWHERE BECAUSE IT RUNS ON A DUAL POWER SUPPLY, AND YOU CAN USE IT IN CARS ALSO, THE MACHINE RUNS ON BOTH DC12 V AND AC 220 V @ 50 HZ

INBUILT BATTERY OPTION:

THE SUCTION MACHINE HAS INTEGRATED 12 V LITHIUM BATTERIES, THE BATTERIES CAN RUN FOR UP TO 60 MINUTES ON A SINGLE FULL CHARGE, WHICH HELPS TO RUN THE MACHINE SEAMLESSLY IN AN EMERGENCY.

HIGH-PERFORMANCE SUCTION:

THE MACHINE PROVIDES A STRONG 25 LPM (LITERS PER MINUTE) SUCTION CAPACITY, WITH AN OIL-FREE, MAINTENANCE-FREE SUCTION PUMP TO REMOVE MUCUS AND PHLEGM.

QUIET OPERATION:

THE 7E-G1 GUARANTEES TO MAINTAIN A CALM ATMOSPHERE IN YOUR HOSPITAL, WHICH IS HELPFUL FOR THE PATIENTS, WITH AN OPERATING NOISE LEVEL OF LESS THAN 50 DB

ADJUSTABLE PRESSURE:

TO MAINTAIN THE NEGATIVE PRESSURE UP TO 80KPA, A REGULATOR IS PROVIDED ON THE FRONT PANEL -

ROBUST DESIGN:

7E-G1 PORTABLE SUCTION MACHINE HAS AN ABS BODY THAT IS LIGHTWEIGHT AND STRONG, MAKING IT MORE PORTABLE AND LONG-LASTING

ENHANCED SAFETY:

FOR A COMPLETE HYGIENE AND SAFETY MACHINE COMES WITH A COMPLETE SETUP INCLUDING, A TUBING KIT, EXTRA BACTERIAL FILTER, AND 1000 ML POLYCARBONATE JAR WITH OVERFLOW PREVENTION FEATURE.

CONVENIENT USE:

ITS LONG-LIFE BATTERY, DETACHABLE POWER CORD, AND LIGHTWEIGHT BODY MAKE IT MORE PORTABLE AND EASY TO CARRY

Jar




Grab
Handle



CE

Procare

Address: Gate no 197, Dehu Alandi road,
Jadhavwadi, Chikali Pune,
Maharashtra-411062
Mob no: +91 8377832832
Web Site : www.procaredmedfurn.com
GST NO: 27ABCFP8176F1ZY

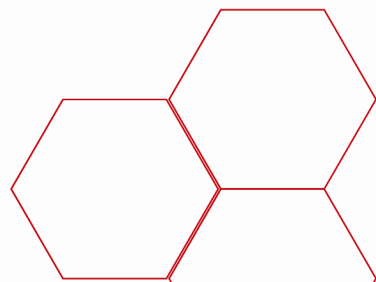

Designs and specifications are subject to change without notice. All images of the product, its components and accessories may not reflect the final product. The images may reflect optional equipment and are subject to change without notice. All information presented herein is based on current data or design goals available at the time of posting and is subject to change without notice. Some equipment and/or features may be unavailable. Actual results of performance and other specifications may differ or vary with production models and may depend on selected models and/or accessories.

Procare